

TAZ 19

Benchtype Fume Adsorber



- *Compact all-round solution for professional labs*
- *Safe and hygienic working place*
- *With down draught extraction method effectively and efficiently protecting user's health*
- *Special hood provides absolute safety against contamination of dangerous fumes*

TAZ 19

Benchtype Fume Adsorber

Safe working place with down-draught extraction and indication of filter saturation by means of a working hours counter for optimum user protection against odor and fume contamination. The large working area allows a wide variety of uses, viz. cutting-up specimens, mounting sections from xylene, and the preparation of reagents. Compared to other instruments available on the market, the **TAZ 19** protects the user from the contamination of dangerous fumes before they pass the activated carbon filter. The suction output can easily be adjusted steplessly according to the requirements of the individual laboratory.

The Benchtype Fume Adsorber **TAZ 19** is offered as three different models:

1. with filter casing including one package of charcoal for xylol
2. same as model 1, adapted for the use of formalin
3. with connection to an external exhaust facility

As an option, a special hood facilitating the exhaust of solvent fumes can be ordered separately (see image below).


Technical Data:

Dimensions :

Overall (W/D/H):	485 x 620 x 405 mm
Base overall (W/D):	485 x 620 mm
Working plate (W/D/H):	310 x 590 x 40 mm
Weight:	20 kg
Power supply:	115-230 V / 50/60 Hz / 160 VA

TAZ 19.000 with special hood

Improved adsorption of solvent fumes

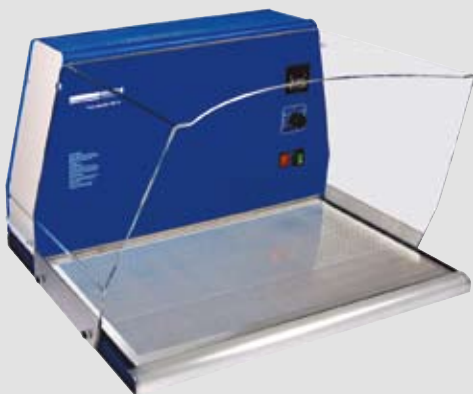

TAZ 19.400 Exhaust (External Connection)

Instrument without charcoal filter and filter casing for access to existing exhaust

07-2008 / Subject to change without prior notice.

MEDITE 

MEDITE GmbH

Wollenweberstr. 12

D-31303 Burgdorf/Germany

Tel. +49 (0) 51 36/88 84-0

Fax +49 (0) 51 36/88 84-55

www.medite-group.com

Distributor:

Made in Germany